


EzWay™ Protein-PreBlue Ladder, 11 bands

- 1. Catalog No.** K18010
- 2. Quantity** 500 ul (50 mini-gels, ready-to-use)
- 3. Storage** Store at 4°C (Store at -20°C for longer storage)
- 4. Description** EzWay™ Protein-PreBlue Ladder is a mixture of highly purified recombinant proteins coupled to a blue dye that resolves to 11 sharp bands during electrophoresis. It allows easy visualization of the MW of proteins and evaluation of the blotting efficiency.
- Pre-stained bands for evaluation of the transfer efficiency
 - Contains 11 size protein bands
 - Extra intense on 2-bands for easy identification (30 kDa, 70 kDa)
 - No lot-to-lot variation
 - Ready-to-use : no need to boil.

- 5. MW Range** 10 - 210 kDa
- 6. No. of Applications** 50 mini-gels (10 ul each)

7. Contents

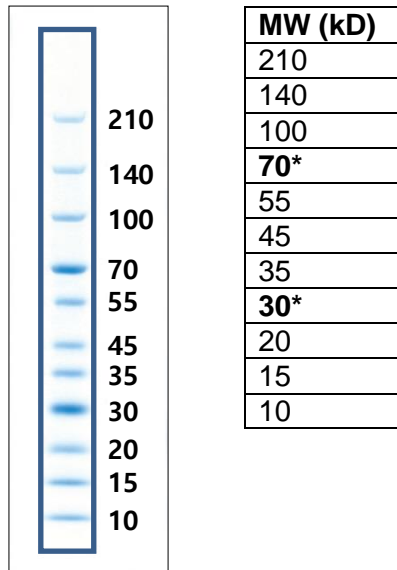

- 8. How to use** Before using, dissolve sediment warming by hand for about 1 minute and mixing lightly, and then load directly 10 ul on to the mini-gel.

- 9. Note** This marker is not suitable for determination of exact molecular weight. For accurate molecular weight estimation of proteins, use unstained Protein marker.(ex/ K08000 EzWay™ Protein-Blue MW Marker)

This marker is not suitable for Aspartate buffer system